

LipoFilter®

Fat Processing System



MICROAIRE®
For Surgery. For Life.™

Gravity-Decant Fat Processing System

LipoFilter® Fat Processing System

LipoFilter is a sterile-packaged, large-volume, fat processing system.

Designed for in-line compatibility with the MicroAire PAL® System*, as well as any traditional suction-assisted liposuction systems (SAL), the LipoFilter offers an innovative approach to filtering adipose tissue using a trouble-free, disposable canister employing a gravity-filtration system.

With a usage capacity of 2500mL, the LipoFilter can meet all of your large volume harvesting needs at an affordable price.

A removable cap on the top of the canister allows easy filling with prewash fluid prior to harvesting, and a unique vent tube provides complete evacuation of fluids after separation is achieved. The vent system also protects fat cells from unnecessary G-forces that occur with other methods.

*Learn more about the PAL System and view the brochure online on www.microaire.com.



LipoFilter is a closed, sterile, and straightforward system that meets all of your large volume harvesting needs.



LipoFilter Aspiration Canister

ASP-CAN-2S (Includes 60cc Toomey Syringe)

ASP-CAN-NT (Syringe Not Included)

The straightforward design of the LipoFilter provides an effective, affordable fat processing system focused on ease of setup and ease of use.

- 3-liter (3000ml) canister
- Compatible with any aspiration system
- Includes filtering canister tubing
- Closed, sterile system



LipoFilter Stand
ASP-CAN-2R


LipoFilter Clamp
ASP-CAN-2C



Adapter, Toomey-to-Luer
Syringe | ASP-ADP-2



Adapter, Toomey-to-Luer
Needle | ASP-ADP-3

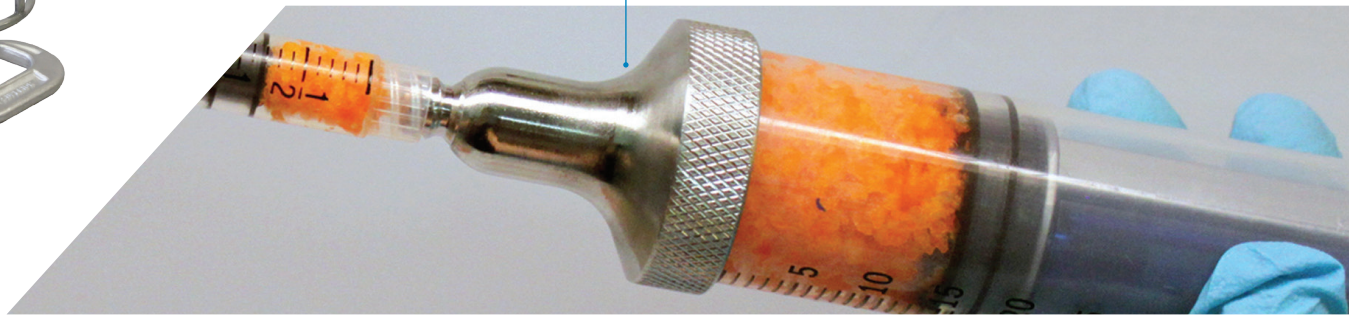

The LipoFilter is designed to provide true ease of use for the surgeon.



NOTE: Syringe not included with ASP-CAN-NT

Ordering Information

LipoFilter® Part Numbers

ASP-CAN-2S	LipoFilter - Aspiration Canister 3000ml, sterile (includes filtering canister tubing and 60cc syringe)
ASP-CAN-NT	LipoFilter - Aspiration Canister 3000ml, sterile (includes filtering canister tubing, syringe not included)
ASP-CAN-2R	LipoFilter Stand
ASP-CAN-2C	LipoFilter Clamp
ASP-ADP-2	Adapter - Toomey-to-Luer Syringe
ASP-ADP-3	Adapter - Toomey-to-Luer Needle

Forhandler: Medistim Norge - epost: norge@medistim.com

Not all fat filters are created equal.

The LipoFilter is designed and manufactured by MicroAire to meet strict FDA standards, because we have your patient's safety in mind.

LipoFilter

FDA 510(K) Approval



FDA-Registered Manufacturing



ISO-Certified Facility



Validated Vendor Supply Chain



Assembled in Class-10,000 Clean Room



ISO 11137-2 Sterilization Validation



Made in the USA



Learn more about our products, including the PAL® System, on www.microaire.com.



MicroAire Surgical Instruments
3590 Grand Forks Boulevard
Charlottesville, VA 22911 USA
800 722 0822 or +434 975 8000
www.microaire.com

© 2022 MicroAire Surgical Instruments, LLC
All rights reserved. LIT-ASP-CAN Rev D



A COLSON ASSOCIATE